



EPITOME OF PERFECTION

Bottle Top Dispenser

(8182)63-90-72
(7172)727-132
(8512)99-46-04
(3852)73-04-60
(4722)40-23-64
(4832)59-03-52
(423)249-28-31
(844)278-03-48
(8172)26-41-59
(473)204-51-73
(343)384-55-89
(4932)77-34-06

(3412)26-03-58
(395)279-98-46
(843)206-01-48
(4012)72-03-81
(4842)92-23-67
(3842)65-04-62
(8332)68-02-04
(861)203-40-90
(391)204-63-61
(4712)77-13-04
(4742)52-20-81

(3519)55-03-13
(495)268-04-70
(8152)59-64-93
(8552)20-53-41
(831)429-08-12
(3843)20-46-81
(383)227-86-73
(3812)21-46-40
(4862)44-53-42
(3532)37-68-04
(8412)22-31-16

(342)205-81-47
- - (863)308-18-15
(4912)46-61-64
(846)206-03-16
- (812)309-46-40
(845)249-38-78
(8692)22-31-93
(3652)67-13-56
(4812)29-41-54
(862)225-72-31
(8652)20-65-13

(3462)77-98-35
(4822)63-31-35
(3822)98-41-53
(4872)74-02-29
(3452)66-21-18
(8422)24-23-59
(347)229-48-12
(4212)92-98-04
(351)202-03-61
(8202)49-02-64
(4852)69-52-93

(996)312-96-26-47

(495)268-04-70

(772)734-952-31

<https://aculab.nt-rt.ru/> || axf@nt-rt.ru

Product Features

Satisfaction Guaranteed

Experience how it feels to be completely satisfied with a device which is designed and built to offer unmatched features that make using it a complete pleasure. Made from the absolute best quality raw materials, Bottle Top Dispensers promise years of durability while offering unwavering precision.



Easy To Dismantle,
Clean And Assemble
The Lower Assembly
With Upgraded Design




Volume Lock
Ensures Stopping Of Any
Accidental Discharge



Unique Piston Mechanism
Allows Cleaning And Servicing
Of Piston Cylinder Without
Disturbing Calibration



The Device Comes With
4 Additional Adapters For
Most Common Bottle Sizes
For Enhanced Convenience

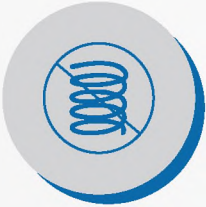

Springless Valve Design
Offers High-chemical
Resistance & A Longer-life
With Minimal Maintenance



Versatile Telescopic Filling
Tube Compatible With Different
Sizes Of Reagent Bottles From
170 MM To 300 MM



Discharge Tubes With
Recirculation Valve Eliminates
Sample Loss & Saves Precious
Reagents & Samples

Convenience at your fingertips

A Superior Way To Work

While ensuring impeccable accuracy, Bottle Top Dispenser is also designed to offer a superior user-experience through various features which make the device extremely easy and convenient to use. Thus, the secret behind its wild popularity is its immense user-friendliness made possible by its highly effective features.



Discharge Tube With 360° Rotation, Makes Dispensing In Multiple Vessels Convenient



Dual Graduation Marks Enable Precise Volume Setting & Clear Visibility From Both Sides.



The Device Is Easy To Recalibrate With A Single Tool.



Calibrated As Per ISO 8655 Standard In A Lab Accredited To ISO 17025 Standard. Each Dispenser Come With Individual Calibration Report.

Ergonomic design with smooth movement for convenient sampling & low pressure dispensing



Snap-fit Mounting Of Discharge-tube For Uncomplicated Mounting



Safety Locking Feature Prevents Accidental Flow Of Reagents



Ordering Information

Bottle Top Dispenser Technica | Specifications

Catalog No. With Re-circulation	Catalog No. Without Re-circulation	Volume Range ml	Increment ml	Inaccuracy (±) %	Inaccuracy (±) ml	Imprecision (±) %	Imprecision (±) ml
ADR - 2.5	AD - 2.5	0.25 to 2.5	0.05	0.6	0.015	0.2	0.005
ADR - 05	AD - 05	0.5 to 5	0.1	0.5	0.025	0.2	0.01
ADR - 10	AD - 10	1 to 10	0.2	0.5	0.05	0.2	0.02
ADR - 25	AD - 25	2.5 to 25	0.5	0.5	0.125	0.2	0.05
ADR - 30	AD - 30	2.5 to 30	0.5	0.5	0.15	0.2	0.06
ADR - 50	AD - 50	5 to 50	1.0	0.5	0.25	0.2	0.1
ADR - 60	AD - 60	5 to 60	1.0	0.5	0.3	0.2	0.12
ADR - 100	AD - 100	10 to 100	2.0	0.5	0.5	0.2	0.2

Accessories

Description	Cat#	Color
Bottle Cap 28mm	ADA-01	
Bottle Cap 38mm	ADA-02	
Bottle Cap 40mm	ADA-03	
Bottle Cap 45mm	ADA-04	
Calibration Tool	ADA-05	
Recirculation Tube	ADA-06	
Telescopic Fill Tube 2.5-60ml	ADA-07	

(8182)63-90-72
 (7172)727-132
 (8512)99-46-04
 (3852)73-04-60
 (4722)40-23-64
 (4832)59-03-52
 (423)249-28-31
 (844)278-03-48
 (8172)26-41-59
 (473)204-51-73
 (343)384-55-89
 (4932)77-34-06

(3412)26-03-58
 (395)279-98-46
 (843)206-01-48
 (4012)72-03-81
 (4842)92-23-67
 (3842)65-04-62
 (8332)68-02-04
 (861)203-40-90
 (391)204-63-61
 (4712)77-13-04
 (4742)52-20-81

(996)312-96-26-47

(3519)55-03-13
 (495)268-04-70
 (8152)59-64-93
 (8552)20-53-41
 (831)429-08-12
 (3843)20-46-81
 (383)227-86-73
 (3812)21-46-40
 (4862)44-53-42
 (3532)37-68-04
 (8412)22-31-16

(495)268-04-70

(342)205-81-47
 (863)308-18-15
 (4912)46-61-64
 (846)206-03-16
 (812)309-46-40
 (845)249-38-78
 (8692)22-31-93
 (3652)67-13-56
 (4812)29-41-54
 (862)225-72-31
 (8652)20-65-13

(772)734-952-31

(3462)77-98-35
 (4822)63-31-35
 (3822)98-41-53
 (4872)74-02-29
 (3452)66-21-18
 (8422)24-23-59
 (347)229-48-12
 (4212)92-98-04
 (351)202-03-61
 (8202)49-02-64
 (4852)69-52-93